

# Accuflow<sup>®</sup> Thala Pump



## vTitan Corporation Private Limited

New Woodlands Building, Estancia IT Park,  
Plot No. 140,151 Vallancheri Village,  
Chengalpattu District,  
Tamil Nadu - 603202, India.

[www.vtitan.com](http://www.vtitan.com)  
[sales@vtitan.com](mailto:sales@vtitan.com)

©2023 vTitan Corporation Pvt. Ltd.  
980-1062 Rev 05







# A Helping Hand to Fight Thalassemia...


## Why Thala Pump :

Thalassemia, an inherited blood disorder characterized by less haemoglobin often requires frequent blood transfusions, resulting in iron overload. Thala pump infuses chelating agents to remove the excess iron.

## Product Features:

-  Compact and Easy to use
-  Simple priming procedure
-  Lock functionality
-  Visual and audio alarms

## 3-Step Operation

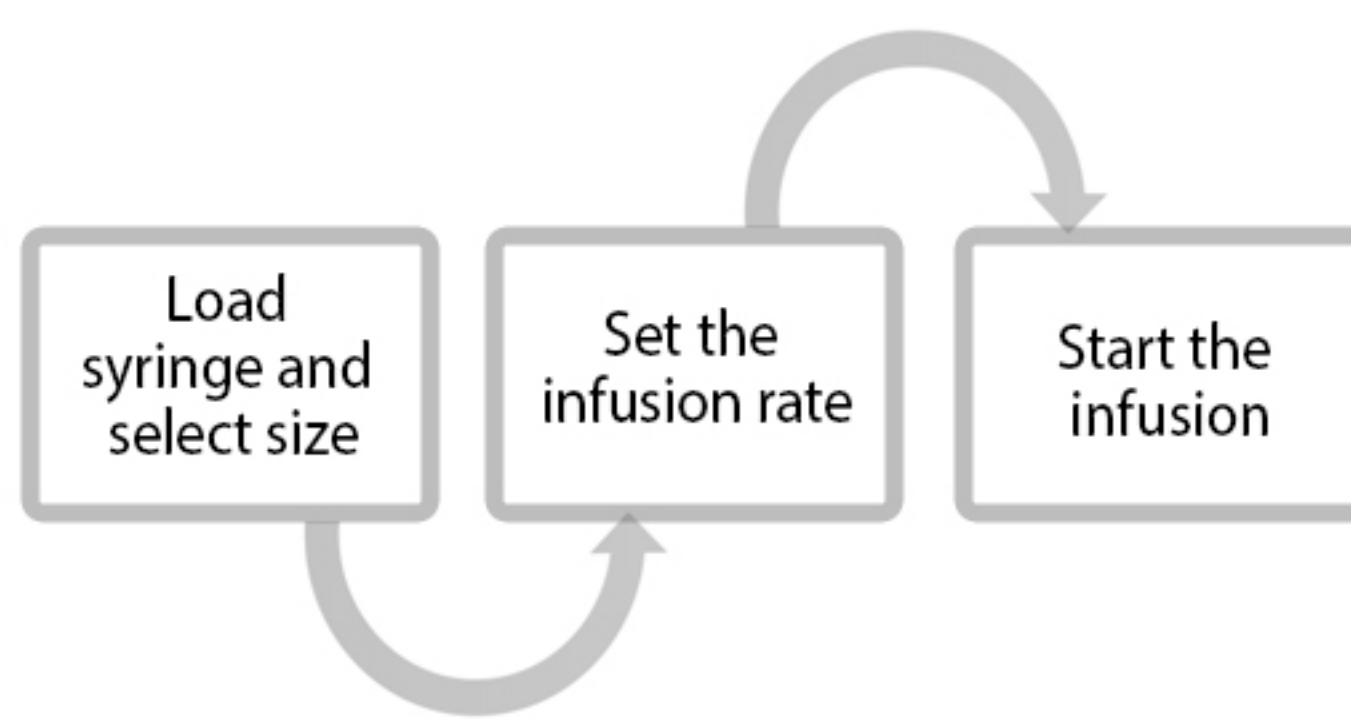

## Product Specifications:

Dimension (in mm)	172 x 36 x 51 (L x W x H)
Weight (in grams)	185 (approx.)
Compatible syringe sizes	5, 10, 20 mL
Display	0.96 inch OLED
Flow rate range	0.1 to 70.0 mL/hr
Accuracy	± 2%
Priming function	Available
Battery type	9V Li-Ion battery (Rechargeable)
Battery life	40 hrs (approx.)
Alarms	Occlusion, Syringe near empty Syringe empty, Syringe disengaged, Battery low, Battery depleted, Unit malfunction
Accessories	Lock box, pouch, charger and adapter
Warranty	1 year